



# The ALPHA

Dolomite's Alpha walkers offer walking support for people with arthritis, neurological disorders, Parkinsons, stroke survivors or people with limitations in wrist or hand movement. A variety of adjustments can be made for arm supported walking. You can use forearms to propel the walker when there is limited hand strength. Arm height adjusts over 9" and the rotation of the arms adjusts 20°.

## Components

Dolomite uses a welded steel frame, no hardware to loosen at joints. Alpha frame has been redesigned wider and deeper to accommodate the higher center of gravity caused by the arm supports. Alpha has 8" high friction, indoor/outdoor tires. Rear of padded forearm supports is recessed to relieve pressure on the elbow. An Allen key is included for adjustments to the arm supports. Step-less arm height adjustment has measurement marks printed on tube to duplicate settings. Flip-up seat gives the user the ability to "walk inside" the frame for better gait training.

## Braking

Alpha Advanced brakes use a padded bar grip. Brake handles can be adjusted so that little force or leverage is needed. A ball joint on allows you to change the direction of brake operation. Alpha Basic has a loop brake handle.

## Standard Features


- Folding frame
- Flip-up seat for gait training and greater reach
- Curb climber on rear wheel assists when lifting front of walker
- Removable basket with carry handle
- Magnetic security latch prevents accidental closure
- Attached Allen key for arm adjustments



## CLARKE HEALTH CARE PRODUCTS, INC.

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

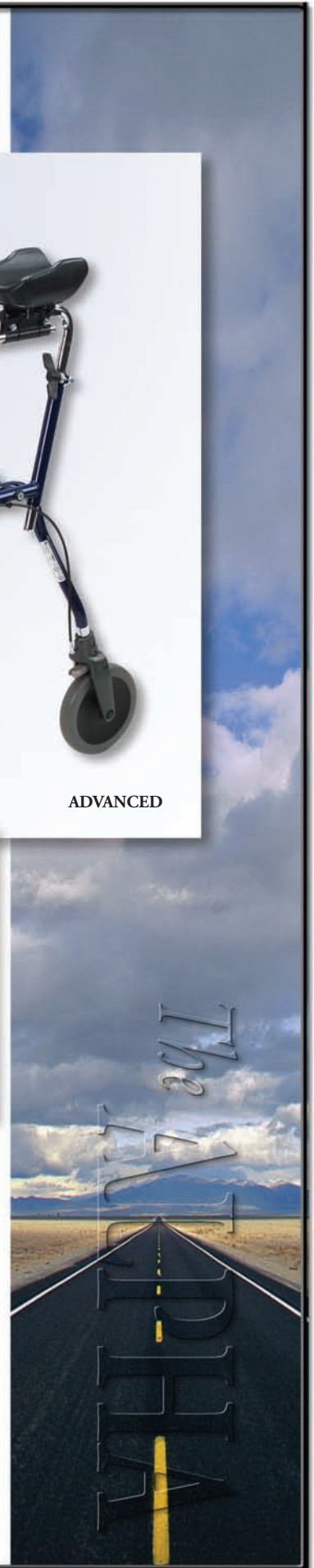



# The ALPHA

## Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



FOOD & BEVERAGE TRAY



CANE HOLDER



PLATFORM SUPPORT

Dolomite has extensive accessory offerings. Accessories can enhance the use of your walker; for instance, the Tray gives you a platform to carry food and the Oxygen Tank Holder will safely support a C or D tank. You can secure a cane or crutch to the side of the walker with the Cane Holder. Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes control standard brakes from one handle. Platform Support replaces the arm troughs; large padded surface to lean onto. Use with either brake system.

SPECIFICATIONS	ALPHA ADVANCED / BASIC
Overall length	27.5"
Overall width	25.5"
Distance between handles	12"-14"
Distance between rear wheels	20"
Seat height	24.5"
Minimum handle height	37.5"
Maximum handle height	47"
Capacity	275 lbs
Walker weight	28 lbs.
Product Number: Advanced	D12220
Product Number: Basic	D12230

### CLARKE HEALTH CARE PRODUCTS, INC.

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

