

AQUATEC®

Ergo

ACCESSORIES AND ADJUSTMENTS



Accessories

Seating and backrests ▼ ▶



A1603142

- ▶ Standard soft seat overlay. With key shaped cut out.



A1603143

- ▶ Soft seat insert. Can be used with or without the soft seat when commode function is not required.




A1603145

- ▶ Small soft seat overlay. With small, oval opening for smaller and young adults.



A1603152

- ▶ Small soft seat insert. Can be used with the small soft seat when the commode function is not required.

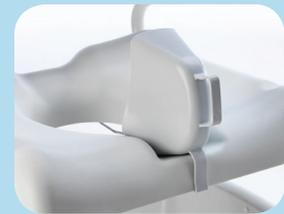

A1603155

- ▶ Splash guard. Clips on easily for use with the standard plastic seat, standard soft seat overlay and small soft seat overlay.



A1531259

- ▶ Ergonomic soft seat overlay, extra soft, textured surface for pelvic stability. 4.6" thick.



A1544469

- ▶ Variable soft seat. Hygiene cut out can be positioned in all directions.



A1603154

- ▶ Universal soft seat overlay. Comfortable closed seat surface.



A1603577

- ▶ XL back. Increases the width between armrests by 3". Included with ErgoXL model, can be retrofitted to standard Ergo model.



A1575299

- ▶ Solid fixed backrest. Offers a stable upright position and is easy to clean.



A1470086

- ▶ Wider XL armrest. Increases distance between armrests by 1.5" each side.



A1558769

- ▶ Soft backrest cushion. For additional stability and comfort. Will also benefit smaller users and young adults. Helps to decrease seat depth.



A1603606, A1604136, A1604016

- ▶ Comfort backrest. For additional comfort, adjustable tension, quick-release buckles. Machine washable.

Positioning ▾



A1535077

- ▶ Side lateral supports. For safe trunk positioning, Width and height adjustable trunk support.



WDSS100P

- ▶ Swing-away side laterals. Adjustable length, swing-away bracket for transfers.



WDACD100

- ▶ Arm trough. Provides additional padded support for arm.



D10167, A1470078

- ▶ Safety Support bar. Additional safety and reassurance, attaches to arm supports. XL size fits ErgoXL.



A1531613

- ▶ Ergonomic headrest holder. Safe and comfortable with various adjustment positions. Pad not included.



A1531608

- ▶ Ergonomic neckrest. Offers exceptional management of the head position.



A1531465

- ▶ Ergonomic headrest. Larger pad for head positioning.



A1470081, A1470082

- ▶ Padded chest belt, seat belt. For safe positioning., Individually adjustable, machine washable.



A1643572, A1643543

- ▶ Calf support with foot support. Padded ELR both height and angle adjustable.



A1643573, A1643574

- ▶ Calf and amputee support. Padded, both height and angle adjustable. For use with the standard plastic seat, standard soft seat overlay and small soft seat overlay.



BB218, BB220

- ▶ Calf strap. Snaps onto footrest tubes, Can prevent feet from kicking behind footrests.



A1470077

- ▶ Arm rest lock. Locks armrest in up or down position.

Toileting ▶



D10230

- ▶ Commode pan with lid. Ergonomic handle allows caretaker to easily remove from pan holder. Scale inside pan.



A1603992, A1601986

- ▶ Front and back anti-tippers. For additional safety in transfers or positioning.

Adjustments

All Ergo models can be adjusted to suit a range of user sizes.

Armrests

- ▶ Ergonomic, swing away armrests can be adjusted to different heights.

Footrests

- ▶ Height adjustable, swing away and detachable footrests
- ▶ No tools required.

Seat Height

- ▶ Seat to ground height can be adjusted between 18.5" - 23.4"
- ▶ No tools required.

Footplates

- ▶ Angle adjustable, flip-up and textured footplates with heel loops
- ▶ No tools required.



Introduction

This useful guide explains how the Aquatec Ergo Family can be customized to suit a wide range of physical conditions by making simple adjustments or by adding some easy to fit, off-the-shelf accessories.

Meet the Ergo Family...



Ergo SPXL ▲

Plus-size self-propelling model

Ergo XL ▲

Plus-size attendant propelled model

Ergo VIP ▲

Tilt-in-space model

Ergo ▲

Attendant propelled model

Ergo Dual VIP ▲

Tilt-in-space model with backrest recline

Ergo SP ▲

Self-propelling model

Compatibility

✓ - Compatible
 X - Not compatible
 S - Included

Most accessories are cross compatible with all Ergo models.

	Product Code	Ergo	Ergo SP	Ergo XL	Ergo SPXL	Ergo VIP	Ergo Dual VIP
Standard soft seat	A1603142	✓	✓	✓	✓	✓	✓
Soft seat insert	A1603143	✓	✓	✓	✓	✓	✓
Small soft seat	A1603145	✓	✓	✓	✓	✓	✓
Small soft seat insert	A1603152	✓	✓	✓	✓	✓	✓
Splash guard	A1603155	✓	✓	✓	✓	✓	✓
Universal Soft seat	A1603154	✓	✓	✓	✓	✓	✓
Ergonomic Special soft seat	A1531259	✓	✓	✓	✓	✓	✓
Variable soft seat	A1544469	✓	✓	✓	✓	✓	✓
Safety bar	D10167	✓	✓	X	X	✓	✓
Safety bar XL	A1470078	X	X	✓	✓	X	X
Chest belt	A1470081	✓	✓	✓	✓	✓	✓
Pelvic/seat belt	A1470082	✓	✓	✓	✓	✓	✓
Armrest lock	A1470077	✓	✓	✓	✓	✓	✓
Arm Trough	WDACD100	✓	✓	✓	✓	✓	✓
XL Back	A1603577	✓	✓	S	S	✓	X
Fixed/solid backrest	A1575299	✓	✓	X	X	✓	✓
Soft backrest cushion	A1558769	✓	✓	✓	✓	✓	✓
Wider XL armrest	A1470086	✓	✓	X	X	✓	✓
24" wheels (set)	A16379	✓	S	✓	S	✓	✓
Calf and Amputee support right	A1556893	✓	✓	✓	✓	✓	✓
Calf and Amputee support left	A1556894	✓	✓	✓	✓	✓	✓
Calf support w/ foot support right ELR	A1546439	✓	✓	✓	✓	✓	✓
Calf support w/ foot support left ELR	A1546705	✓	✓	✓	✓	✓	✓
Calf strap 18"	BB218	✓	✓	✓	✓	✓	✓
Calf strap 20"	BB220	✓	✓	✓	✓	✓	✓
Side lateral	A1535077	✓	✓	✓	✓	✓	✓
Swing-away side lateral	WDSS100P	✓	✓	✓	✓	✓	✓
Ergonomic headrest holder	A1531613	X	X	X	X	✓	✓
Ergonomic headrest	A1531465	X	X	X	X	✓	✓
Ergonomic neckrest	A1531608	X	X	X	X	✓	✓
Commode pan with lid	D10230	✓	✓	✓	✓	✓	✓
Guide rail for commode pan	A1602078	✓	✓	✓	✓	✓	✓
Front and back anti-tipper	A1603992	X	X	X	X	✓	✓
Front anti-tipper	A1601986	✓	✓	✓	✓	X	X
Comfort soft backrest ErgoVIP	A1603606	X	X	X	X	✓	✓
Comfort soft backrest Ergo	A1604136	✓	✓	X	X	X	X
Comfort soft backrest ErgoXL	A1604016	X	X	✓	✓	X	X

Clarke Healthcare Products
 7830 Steubenville Pike, Oakdale, PA 15071
www.clarkehealthcare.com 888-347-4537

