

**AQUATEC®**

# Ergo VIP

# Ergo Dual VIP



**Modular personal care solutions  
for complex needs**

ERGONOMIC DESIGN | SAFETY | INDEPENDENCE



# Modular personal care solutions for those with more complex seating and positioning needs

## Highly customizable


The **Aquatec Ergo VIP** and **Ergo Dual VIP** are modular tilt in space shower chair commodes designed to suit a wide range of individuals and their postural needs. Both models benefit from a -5° to 40° seat tilt angle, whilst the Dual VIP has an additional backrest recline of 0° to 35° for those with more complex seating and positioning requirements.

As with all Ergo models, the VIP and Dual VIP can be tailored to each individual by making some simple adjustments or by adding some easy to fit, off the shelf accessories. Furthermore, their ergonomic design ensures ease of use for the carer whilst providing optimum safety and comfort to the user.

See full range of accessories online at [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

### NEW Contoured, ergonomic seat

- ▶ Based on the imprints people leave behind when seated to ensure a more upright and stable seating position.
- ▶ A 5° seat tilt to lift the knees into the natural squat position for toileting.
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position for toileting and showering.

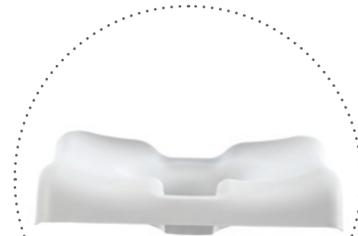

### NEW Ergonomic tilt

- ▶ The new pivot and slide mechanism means that during tilt, the chair now moves downwards and forwards which keeps the user's center of gravity close to that of the chair. This makes the tilt operation lighter and more intuitive for the carer, creates a smoother and more comfortable experience for the user and has reduced the overall footprint of the product without affecting its stability.



**NEW**

FEATURES

**-5° to 40°**

**1** **Extremely stable**

- ▶ Enhanced stability for additional safety and reassurance

**2** **Ergonomic tilt**

- ▶ Ergonomic tilt with weight shifting abilities

**3** **Innovative ergonomic seat**

- ▶ Contoured, ergonomic seat plate encourages a more upright seating position for enhanced independence and comfort

**4** **Smaller footprint**

- ▶ Optimized wheel base makes it easier to maneuver in small spaces



## Other features



▲ Stainless steel corrosion free frame.



▲ Extremely stable even at maximum tilt.



▲ Machine washable tension adjustable backrest.



▲ Armrest height, footrest height and seat height can be easily adjusted.



▲ 5° - 35° Backrest Recline  
**Ergo Dual VIP Only**



# Meet the Ergo VIP and Ergo Dual VIP ...

---

Modular tilt in space models to suit a wide range of individuals and their postural needs.



***Ergo VIP***

▲ -5°-40° seat tilt angle

---



***Ergo DualVIP***

▲ -5°-40° seat tilt angle with an additional  
backrest recline of 5° - 35°

---

## Technical Data ▼

For more comprehensive information about this product, including the product's user manual, please see your Clarke Health Care rep or Aquatec dealer.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

<b>Ergo VIP</b>	18.7"	18.7"	22"	39"	17.75"
<b>Ergo Dual VIP</b>	18.7"	18.7"	22"	39"	17.75"



Tilt angle seat ▼



Tilt angle back ▼



Weight ▼



Load capacity ▼

<b>Ergo VIP</b>	-5° - 40°	-	51.8 lb.	330 lb.
<b>Ergo Dual VIP</b>	-5° - 40°	5° - 35°	62.4 lb.	330 lb.

## Frame color ▼



White

## Regulatory data ▼




**Clarke Healthcare Products**  
 7830 Steubenville Pike, Oakdale, PA 15071  
[www.clarkehealthcare.com](http://www.clarkehealthcare.com) **888-347-4537**



AQUATEC®

# Ergo Family

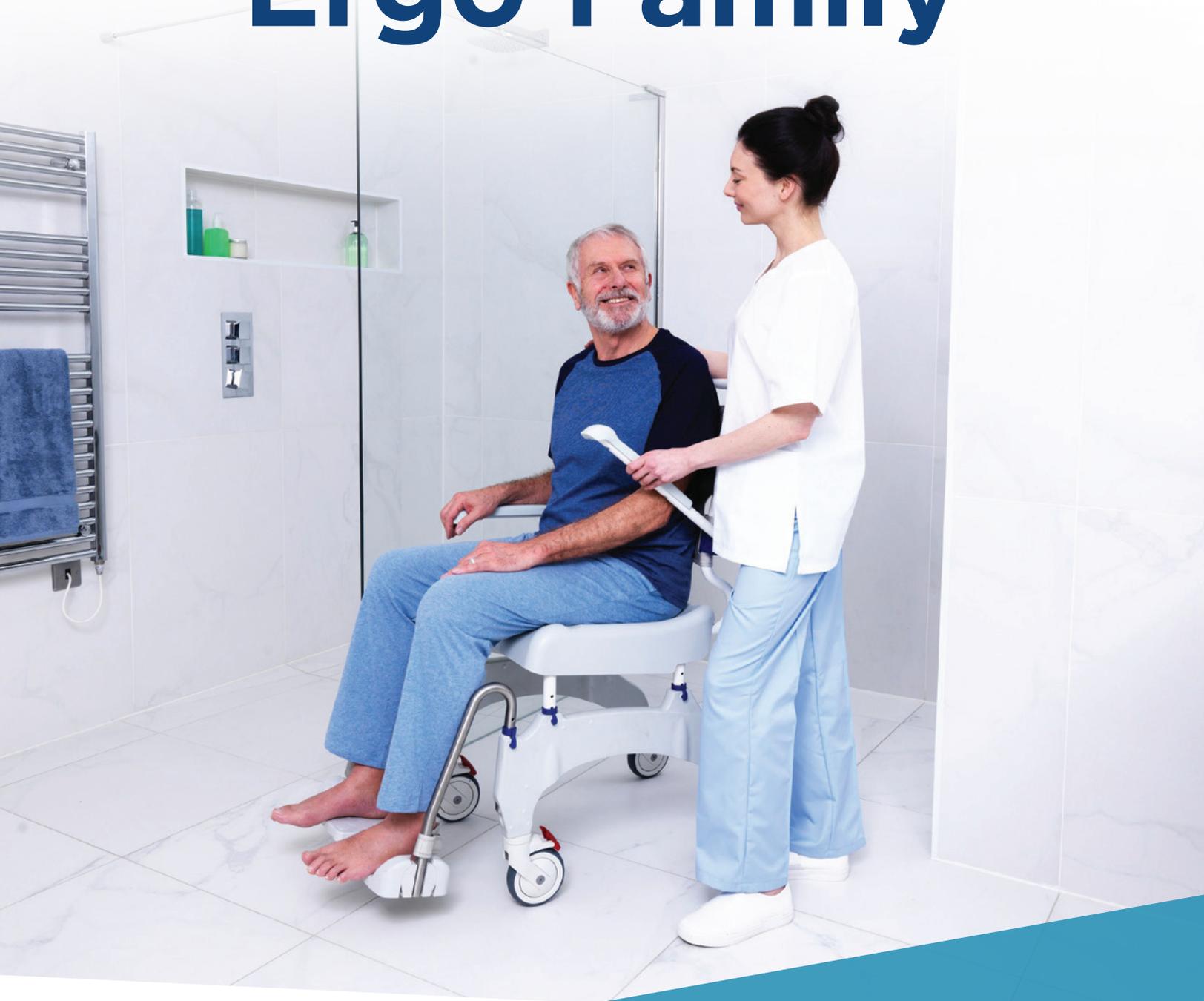

**Modular personal care solutions**

SAFE | INDEPENDENT | ERGONOMIC



# Customer insight-led design to enhance independence, safety and flexibility.

Personal care is highly private activity and is an important part of a person's life and wellbeing. The **Aquatec Ergo Family** of shower chair commodes have been designed with user comfort, independence and dignity in mind. Next generation of the well-established Aquatec shower chairs, the Ergo family comes packed with a range of new features, designed to make personal care activities easier for both the user and carer.

### **Safety, independence and dignity**

A newly contoured seat surface has been developed to encourage a more upright and stable seating position, promoting user independence and overall

comfort. Enhanced stability at the front of the base allows for a smaller footprint and offers user added safety and reassurance when seated in the chair. Furthermore, the optimized wheelbase makes it even easier for a carer to maneuver the chair in small spaces as well as allowing for a better fit over standard toilets for independent toileting.

### **Modular and flexible design**

The modular design of the Ergo family means that the chairs can be adapted to suit a wide range of user needs by making some simple adjustments, or by adding some of our easy to fit, off-the-shelf accessories.

## **OTHER** features include:

- ▶ Stable stainless steel frame
- ▶ Machine washable tension adjustable back with quick-release buckles
- ▶ Easy disassembly for transportation and storage



# NEW FEATURES

## 1 Innovative ergonomic seat

▶ Insight-led design encourages a more upright seating position for enhanced independence and comfort

## 2 Greater adaptability

▶ With an increased height adjustment range on the footrests and armrests

## 4 Compact and easy to handle

▶ Optimized wheel base makes it easy to maneuver in small spaces


## 3 Extremely stable

▶ Enhanced stability at the front for additional safety and reassurance



## NEW Contoured, ergonomic seat plate

- ▶ Based on the imprints people leave behind when seated
- ▶ A 5° seat tilt lifts the knees into a natural squat position for toileting
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position

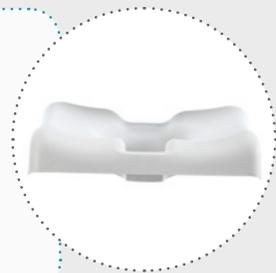

## Highly Configurable

Available with a extensive range of positioning, safety, toileting and seating accessories that can be used to support a wide range of conditions and postural needs.


See full range of accessories online at [www.clarkehealthcare.com](http://www.clarkehealthcare.com)



Positioning



Safety

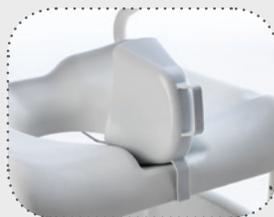

Toileting



Seating and Backrests

## Tilt in space models:

### Ergo VIP and Ergo Dual VIP

- ▶ Both have a -5-40° seat tilt angle with a smart weight shift mechanism that ensures optimal stability on a small footprint, whilst the Ergo Dual VIP has an additional backrest recline of 5°-35°



# Introducing the **Ergo** Family

From a standard transit model, through to tilt and recline, meet the flexible range of shower chairs that offer personal care solutions for a wide range of physical conditions.



## **Ergo SPXL** ▲

An entry level plus-size model with a maximum user weight of 400 lb.

## **Ergo XL** ▲

An entry level plus-size attendant propelled model with a maximum user weight of 400 lb.

## **Ergo VIP** ▲

Designed for individuals who require slightly more assistance with postural control when seated

## **Ergo** ▲

An entry level attendant propelled model

## **Ergo Dual VIP** ▲

Designed for individuals with complex seating and positioning needs

## **Ergo SP** ▲

An entry level self-propelling model

## Technical Data ▼

For more comprehensive information about this product, including the product's user manual, please visit [www.clarkehealthcare.com](http://www.clarkehealthcare.com).



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

	Seat width ▼	Seat depth ▼	Total width ▼	Total depth ▼	Width between armrests ▼
<b>Ergo</b>	18.7"	18.7"	22"	35.3"	17.75"
<b>Ergo XL</b>	18.7"	18.7"	25"	35.3"	20.87"
<b>Ergo SP</b>	18.7"	18.7"	26"	41"	17.75"
<b>Ergo SPXL</b>	18.7"	18.7"	26"	41"	20.87"
<b>Ergo VIP</b>	18.7"	18.7"	22"	39"	17.75"
<b>Ergo Dual VIP</b>	18.7"	18.7"	22"	39"	17.75"



Tilt angle seat ▼



Tilt angle back ▼



Weight ▼



Load capacity ▼

	Tilt angle seat ▼	Tilt angle back ▼	Weight ▼	Load capacity ▼
<b>Ergo</b>	-	-	37.5 lb.	330 lb.
<b>Ergo XL</b>	-	-	39.7 lb.	400 lb.
<b>Ergo SP</b>	-	-	50.7 lb.	330 lb.
<b>Ergo SPXL</b>	-	-	52.9 lb.	400 lb.
<b>Ergo VIP</b>	-5° - 40°	-	51.8 lb.	330 lb.
<b>Ergo Dual VIP</b>	-5° - 40°	5° - 35°	62.4 lb.	330 lb.

## Frame color ▼



White

## Regulatory data ▼



**Clarke Healthcare Products**  
 7830 Steubenville Pike, Oakdale, PA 15071  
[www.clarkehealthcare.com](http://www.clarkehealthcare.com) **888-347-4537**



**AQUATEC®**

# Ergo

**ACCESSORIES AND ADJUSTMENTS**



# Accessories

## Seating and backrests ▼ ▶



**A1603142**

- ▶ Standard soft seat overlay. With key shaped cut out.



**A1603143**

- ▶ Soft seat insert. Can be used with or without the soft seat when commode function is not required.



**A1603145**

- ▶ Small soft seat overlay. With small, oval opening for smaller and young adults.



**A1603152**

- ▶ Small soft seat insert. Can be used with the small soft seat when the commode function is not required.



**A1603155**

- ▶ Splash guard. Clips on easily for use with the standard plastic seat, standard soft seat overlay and small soft seat overlay.



**A1531259**

- ▶ Ergonomic soft seat overlay, extra soft, textured surface for pelvic stability. 4.6" thick.



**A1544469**

- ▶ Variable soft seat. Hygiene cut out can be positioned in all directions.



**A1603154**

- ▶ Universal soft seat overlay. Comfortable closed seat surface.



**A1603577**

- ▶ XL back. Increases the width between armrests by 3". Included with ErgoXL model, can be retrofitted to standard Ergo model.



**A1575299**

- ▶ Solid fixed backrest. Offers a stable upright position and is easy to clean.



**A1470086**

- ▶ Wider XL armrest. Increases distance between armrests by 1.5" each side.



**A1558769**

- ▶ Soft backrest cushion. For additional stability and comfort. Will also benefit smaller users and young adults. Helps to decrease seat depth.



**A1603606, A1604136, A1604016**

- ▶ Comfort backrest. For additional comfort, adjustable tension, quick-release buckles. Machine washable.

## Positioning ▾



**A1535077**

- ▶ Side lateral supports. For safe trunk positioning, Width and height adjustable trunk support.



**WDSS100P**

- ▶ Swing-away side laterals. Adjustable length, swing-away bracket for transfers.



**WDACD100**

- ▶ Arm trough. Provides additional padded support for arm.



**D10167, A1470078**

- ▶ Safety Support bar. Additional safety and reassurance, attaches to arm supports. XL size fits ErgoXL.



**A1531613**

- ▶ Ergonomic headrest holder. Safe and comfortable with various adjustment positions. Pad not included.



**A1531608**

- ▶ Ergonomic neckrest. Offers exceptional management of the head position.



**A1531465**

- ▶ Ergonomic headrest. Larger pad for head positioning.



**A1470081, A1470082**

- ▶ Padded chest belt, seat belt. For safe positioning., Individually adjustable, machine washable.



**A1643572, A1643543**

- ▶ Calf support with foot support. Padded ELR both height and angle adjustable.



**A1643573, A1643574**

- ▶ Calf and amputee support. Padded, both height and angle adjustable. For use with the standard plastic seat, standard soft seat overlay and small soft seat overlay.



**BB218, BB220**

- ▶ Calf strap. Snaps onto footrest tubes, Can prevent feet from kicking behind footrests.



**A1470077**

- ▶ Arm rest lock. Locks armrest in up or down position.

## Toileting ▶



**D10230**

- ▶ Commode pan with lid. Ergonomic handle allows caretaker to easily remove from pan holder. Scale inside pan.



**A1603992, A1601986**

- ▶ Front and back anti-tippers. For additional safety in transfers or positioning.

# Adjustments

All Ergo models can be adjusted to suit a range of user sizes.

## Armrests

- ▶ Ergonomic, swing away armrests can be adjusted to different heights.

## Footrests

- ▶ Height adjustable, swing away and detachable footrests
- ▶ No tools required.

## Seat Height

- ▶ Seat to ground height can be adjusted between 18.5" - 23.4"
- ▶ No tools required.

## Footplates

- ▶ Angle adjustable, flip-up and textured footplates with heel loops
- ▶ No tools required.



# Introduction

This useful guide explains how the Aquatec Ergo Family can be customized to suit a wide range of physical conditions by making simple adjustments or by adding some easy to fit, off-the-shelf accessories.

## Meet the Ergo Family...



### **Ergo SPXL ▲**

Plus-size self-propelling model

### **Ergo XL ▲**

Plus-size attendant propelled model

### **Ergo VIP ▲**

Tilt-in-space model

### **Ergo ▲**

Attendant propelled model

### **Ergo Dual VIP ▲**

Tilt-in-space model with backrest recline

### **Ergo SP ▲**

Self-propelling model

# Compatibility

✓ - Compatible  
 X - Not compatible  
 S - Included

Most accessories are cross compatible with all Ergo models.

	Product Code	Ergo	Ergo SP	Ergo XL	Ergo SPXL	Ergo VIP	Ergo Dual VIP
Standard soft seat	A1603142	✓	✓	✓	✓	✓	✓
Soft seat insert	A1603143	✓	✓	✓	✓	✓	✓
Small soft seat	A1603145	✓	✓	✓	✓	✓	✓
Small soft seat insert	A1603152	✓	✓	✓	✓	✓	✓
Splash guard	A1603155	✓	✓	✓	✓	✓	✓
Universal Soft seat	A1603154	✓	✓	✓	✓	✓	✓
Ergonomic Special soft seat	A1531259	✓	✓	✓	✓	✓	✓
Variable soft seat	A1544469	✓	✓	✓	✓	✓	✓
Safety bar	D10167	✓	✓	X	X	✓	✓
Safety bar XL	A1470078	X	X	✓	✓	X	X
Chest belt	A1470081	✓	✓	✓	✓	✓	✓
Pelvic/seat belt	A1470082	✓	✓	✓	✓	✓	✓
Armrest lock	A1470077	✓	✓	✓	✓	✓	✓
Arm Trough	WDACD100	✓	✓	✓	✓	✓	✓
XL Back	A1603577	✓	✓	S	S	✓	X
Fixed/solid backrest	A1575299	✓	✓	X	X	✓	✓
Soft backrest cushion	A1558769	✓	✓	✓	✓	✓	✓
Wider XL armrest	A1470086	✓	✓	X	X	✓	✓
24" wheels (set)	A16379	✓	S	✓	S	✓	✓
Calf and Amputee support right	A1556893	✓	✓	✓	✓	✓	✓
Calf and Amputee support left	A1556894	✓	✓	✓	✓	✓	✓
Calf support w/ foot support right ELR	A1546439	✓	✓	✓	✓	✓	✓
Calf support w/ foot support left ELR	A1546705	✓	✓	✓	✓	✓	✓
Calf strap 18"	BB218	✓	✓	✓	✓	✓	✓
Calf strap 20"	BB220	✓	✓	✓	✓	✓	✓
Side lateral	A1535077	✓	✓	✓	✓	✓	✓
Swing-away side lateral	WDSS100P	✓	✓	✓	✓	✓	✓
Ergonomic headrest holder	A1531613	X	X	X	X	✓	✓
Ergonomic headrest	A1531465	X	X	X	X	✓	✓
Ergonomic neckrest	A1531608	X	X	X	X	✓	✓
Commode pan with lid	D10230	✓	✓	✓	✓	✓	✓
Guide rail for commode pan	A1602078	✓	✓	✓	✓	✓	✓
Front and back anti-tipper	A1603992	X	X	X	X	✓	✓
Front anti-tipper	A1601986	✓	✓	✓	✓	X	X
Comfort soft backrest ErgoVIP	A1603606	X	X	X	X	✓	✓
Comfort soft backrest Ergo	A1604136	✓	✓	X	X	X	X
Comfort soft backrest ErgoXL	A1604016	X	X	✓	✓	X	X

## Clarke Healthcare Products

7830 Steubenville Pike, Oakdale, PA 15071

[www.clarkehealthcare.com](http://www.clarkehealthcare.com) 888-347-4537

