



The LEGACY

Dolomite developed the **Legacy** to provide the most standard features to a variety of users for the best value. Over the years **Legacy** has proven to be a reliable aid for mobility. **Legacy's** steel frame is welded for greater durability, no screws to loosen at joints. The framework is powder-coated for corrosion-resistance and comes in Metallic Blue or Red. Three frame heights accommodate various body sizes. With forward facing handles and flip-up seat you walk within the **Legacy** footprint with frame supporting you on all sides. Handles are height-adjustable for best posture while walking. Push the seat down to hear a click alerting you that the walker is locked open. **Legacy** is easily folded for storage or transportation, lift seat handle to release and fold.


Legacy comes with the patented Dolomite DABS braking system, with adjustable tension for varied grip strength. Brake handles have a wide flange and large loops for comfortable hand positioning. When parked and seated, push down the handle to set parking brakes. The braking mechanism is located within the fork for protection from damage. Optional accessory slow-down brakes offer adjustable tension on the tires for more control on hilly terrain or to even out direction of travel for stroke patients.

Standard Features

- 330 lb. capacity
- Flip-up seat for greater reach or gait training
- Ergonomic and secure hand grips
- Flexible and adjustable back strap
- 8" indoor-outdoor solid tires
- Optional basket with carry handle

Warranty

To the original purchaser, Dolomite provides a limited life-time warranty against manufacturer's defects on the frame, five years on the brake handle and one year on parts. Tires, brake cables grips are considered wearable and are excluded from this warranty.

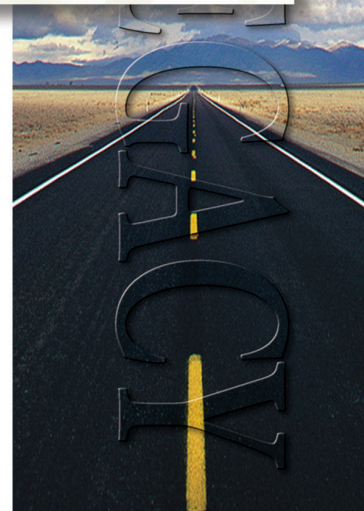

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The LEGACY

Accessories



ANTI-STATIC GRIP



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



CROSS BAR WITH BRAKE



FOOD & BEVERAGE TRAY



CANE HOLDER



TALL HANDLES



UNDER SEAT BAG



SEAT AND BACK PAD

Accessories enhance the use of your rollator, the Tray gives you a surface to carry food, the Oxygen Tank holder will safely support tanks sizes B, C, D, E. A Cane Holder can secure a cane or crutch to the side of our rollator. Slow-Down Brakes have tension control to allow individual setting on each tire. One-Hand Brakes control braking from one handle. Under-Seat bag offers a place to stash items out of sight. Tall Handles add additional height. Cross-Bar combines a forward push bar grip area with one-hand braking.

SPECIFICATIONS	600	520	450
Dist. Between Rear Wheels	19"	19"	19"
Minimum & Maximum Handle Height	33 - 36.5"	31 - 34.5"	28 - 31.5"
Minimum & Maximum Tall Handle Heights	37 - 41"	34 - 38"	32 - 36"
Seat Height from Floor	23.8"	22"	19"
Seat Dimensions	7" Deep x 16.5" Wide		
Overall Width/Length	23"/26"	23"/26"	23"/23"
Width Between Handles	17"	17"	17"
Folded Dimensions	36" Length x 23" Width x 12" Height		
Walker Weight (lbs.)	16	16	16
Weight Capacity (lbs.)	330	330	330
Product Number: Blue	D12050B	D12052B	D12054B
Product Number: Red	D12050R	D12052R	D12054R

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

