



The LEGACY

Dolomite developed the Legacy to accommodate a wide variety of users with the most standard features for the best value. The design incorporates these features: flip-up seat, folding frame, adjustable handles, and wide wheels. With choice of color and an extensive accessory line, the Legacy looks smart strolling down the sidewalk.

Frame

The strongest frame on the market, Dolomite welds the steel frame for greater durability, no screws to loosen at joints. Legacy is available in two colors: Metallic Blue or Candy Apple Red. There are three frame sizes to accommodate individuals of all sizes. Dolomite frames are easily folded for storage or transportation. The handles are height-adjustable in 4 increments.

Braking

Legacy comes with the patented Dolomite Automatic Braking System with adjustable tension for varied grip strength. Brake handles have wide flange and large loops for comfortable hand positioning. When parked and using the seat, push down the handle to set the parking brakes, audible click is heard. The brake mechanism is located within fork for protection from damage. Accessory slow down brake provides adjustable tension for hilly terrain or to even out direction of travel for stroke patients.

Standard Features

Legacy supports 330 pounds.

Other standard features are:

- Flip-up seat for greater reach or gait training
- Seat handle releases to fold walker
- Formed grips for comfortable and secure grip
- Flexible and adjustable back strap
- 8" indoor/outdoor tires

Warranty

To original purchaser, Dolomite provides a lifetime warranty against manufacturer's defects on the frame, 5 years on brake handle and one year on parts. Tires, brake cables and grips are excluded from this warranty.



CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The LEGACY

Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



CROSS BAR WITH BRAKE



FOOD & BEVERAGE TRAY




CANE HOLDER



TALL HANDLES



UNDER SEAT BAG

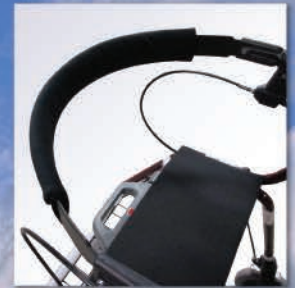

SEAT AND BACK PAD

Accessories can enhance the use of your walker – the Tray gives you a platform to carry food, removable basket with carry handle, the Oxygen Tank Holder will safely support a B, C or D tank. A Cane Holder can secure a cane or crutch to the side of the walker, Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes can control your braking from one handle. Under Seat Bag gives you a place to stash items out of sight. Tall handles add 3” to the handle height. Cross bar combines more hand hold area with one hand braking.

SPECIFICATIONS	SUPER LOW	LOW	STANDARD
Dist. Between Rear Wheels	19"	19"	19"
Minimum & Maximum Handle Height	28 - 31.5"	31 - 34.5"	33 - 36.5"
Minimum & Maximum Tall Handle Heights	31 - 35"	33 - 37"	36 - 40"
Seat Height from Floor	19"	22"	24"
Seat Dimensions	7" Deep x 15" Wide		
Overall Width/Length	23"/23"	23"/26"	23"/26"
Width Between Handles	17"	17"	17"
Folded Dimensions	36" Length x 23" Width x 12" Height		
Walker Weight (lbs.)	16	17	17
Weight Capacity (lbs.)	330	330	330
Product Number: Blue	D12054B	D12052B	D12050B
Product Number: Red	D12054R	D12052R	D12050R

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

