



# The MAXI+

Dolomite designed the Maxi+ for the heavier or larger user. The wider frame provides more seat room and larger walking area but still easily maneuvers in your environment. Over the years Maxi models have been a reliable walking aid. The high-tensile steel frame has welded joints for strength and low maintenance. The powder-coated finish for corrosion-resistance comes in Metallic blue color. Two frame heights offer a choice for best posture. A red safety latch secures the Maxi+ in open position. Lift latch to fold. With forward facing handles and flip-up seat you are walking within the Maxi+ footprint, with frame supporting you on all sides. Ergonomic grips reduce hand fatigue and provide a secure hold.

## Frame


Maxi+ uses patented Dolomite DABS brakes with adjustable tension for varied grip strength. Brake handles have a wide flange and large loops for comfortable hand positioning. The brake mechanism is located within the forks for protection from damage. When parked and seated push down the handle to set the parking brake. Optional accessory slow-down brake offers adjustable tension on the tire for more control on hilly terrain or to even out direction of travel for stroke patients.

## Features

- Extra-wide, 20" between handles
- Large flip-up seat
- 452 lb. capacity
- 8" indoor-outdoor wheels
- Optional basket with carry handle
- Adjustable back support
- Handles adjust over a 7" range

## Warranty

To the original purchaser, Dolomite provides a limited life-time warranty against manufacturer's defects on the frame, five years on brake handle and one year on parts. Tires, brake cables and grips are considered wearable and excluded from this warranty.

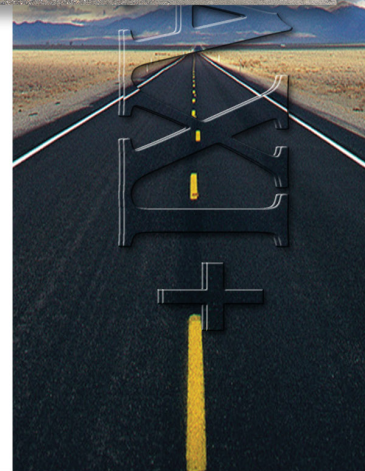

## CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)



# The MAXI+

## Accessories



TRAY



OXYGEN TANK HOLDER



CANE HOLDER



TALL HANDLES  
MAXI+



SLOW DOWN BRAKE



ONE HAND BRAKE



UNDER SEAT BAG



SEAT AND BACK PAD

Accessories can enhance your walker – the Tray can carry food, Oxygen Tank Holder will safely support a B, C, D or E tank. Cane Holder can secure a cane of crutch on the side of the walker. Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes allow braking from one handle. Under Seat Bag gives you a place to stash items out-of-sight. Handle extensions increase height to 41". Large basket for personal items.

### Specifications:

SPECIFICATIONS	550	650
Dist. Between Rear Wheels	22"	22"
Handle Height	29" - 35"	32" - 38"
Tall Handle Height	41"	36" - 42"
Seat Height from Floor	22"	25"
Seat Dimensions	8.5" Deep x 20" Wide	
Overall W x L	27.5" W x 26" L	27.5" W x 29" L
Width Between Handles	20.5"	20.5"
Folded Dimensions	33" Length x 27.5" Width x 11" Height	
Walker Weight (lbs.)	21.6	21
Weight Capacity (lbs.)	452	452
Product Number	D12122	D12120

### CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

